

# NG5000 pureza ultra alta Generador de nitrógenos

Aplicación: Soluciones especializadas


Número de pieza: 10-2310 (230v) - UK | 10-2320 (230v) - EU | 10-2330 (110v) - US



## Description

El NG5000 entrega de nitrógeno de ultra alta pureza en los límites de detección de trazas a través de la tecnología de adsorción por cambio de presión en los flujos de 5 l / min para aplicaciones GC .

- **Tipos de gas** Nitrógeno
- **Compresor integrado:** No

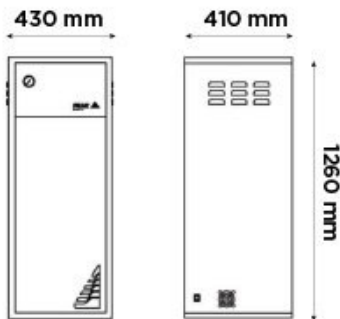

## Key Features

- Adecuado para aplicaciones como GC, DSC, DMA
- Adecuado para aplicaciones que requieren una pureza ultra alta
- Disponible con o sin compresor externo
- Contiene una columna de tamiz molecular de carbón autorregenerante
- El gas se suministra a pedido para que el generador funcione según su horario
- Garantía integral in situ de 12 meses

## Technical Specifications

### NG5000 pureza ultra alta Generador de nitrógenos

Tipo de gas	Nitrógeno
Flujo de gas máximo:	5L/min
Presión máxima de salida:	80psi/5,5bar
Presión máxima de salida:	80psi/5,5bar
Presión mínima de entrada:	120psi/8,3bar
Presión máxima de entrada:	130psi/9bar
Voltaje:	110 / 230 ± 10% VAC
Frecuencia:	50 / 60 Hz
Corriente:	0.5 Amps / 1.2 Amps
Temperatura máxima de funcionamiento:	25°C / 77°F
Acreditaciones:	CE
Tamaño (HxAxD) en mm	1260 x 430 x 410 mm
Tamaño (HxAxD) en pulgadas	49,6 x 16,9 x 16,1 inches
Peso del generador	80kg / 176lbs

# [PEAK Protected]<sup>TM</sup>

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: [www.peakscientific.com/ordering/](http://www.peakscientific.com/ordering/)
- For service plans visit: [www.peakscientific.com/service/service-plans/](http://www.peakscientific.com/service/service-plans/)

**Peak Scientific UK**  
Tel: +44 (0)141 812 8100  
Fax: +44 (0)141 812 8200

**Peak Scientific North America**  
Tel: +1-800-767-6532  
Fax: +1-978-608-9503

**Peak Scientific China**  
Tel: +86 21 5079 1190  
Fax: +86 21 5079 1191  
For a full list of our worldwide office locations, please visit:  
Web: [www.peakscientific.com](http://www.peakscientific.com)  
Email: [discover@peakscientific.com](mailto:discover@peakscientific.com)

Peak Scientific's Quality Management System conforms to: ISO:9001

