

NG3000 pureza ultra alta Generador de nitrógenos

Aplicación: Soluciones especializadas


Número de pieza: 10-2210 (230v) - UK | 10-2211 (230v) - EU | 10-2230 (110v) - US



Description

El NG3000 entrega de nitrógeno de ultra alta pureza en los límites de detección de trazas a través de la tecnología de adsorción por cambio de presión en los flujos de 3 l / min para aplicaciones GC .

- **Tipos de gas** Nitrógeno
- **Compresor integrado:** No

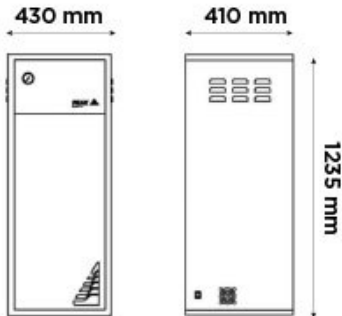

Key Features

- Adecuado para aplicaciones como GC, DSC, DMA
- Adecuado para aplicaciones que requieren una pureza ultra alta
- Disponible con o sin compresor externo
- Contiene una columna de tamiz molecular de carbón autorregenerante
- El gas se suministra a pedido para que el generador funcione según su horario
- Garantía integral in situ de 12 meses

Technical Specifications

NG3000 pureza ultra alta Generador de nitrógenos

Tipo de gas	Nitrógeno
Flujo de gas máximo:	3
Presión máxima de salida:	80psi/5,5bar
Presión máxima de salida:	80psi/5,5bar
Presión mínima de entrada:	120psi/8,3bar
Presión máxima de entrada:	130psi/9bar
Voltaje:	110 / 230 ± 10% VAC
Frecuencia:	50 / 60 Hz
Corriente:	0.5 Amps / 1.2 Amps
Temperatura máxima de funcionamiento:	25°C / 77°F
Acreditaciones:	CE
Tamaño (HxAxD) en mm	1260 x 430 x 410 mm
Tamaño (HxAxD) en pulgadas	49,6 x 16,9 x 16,1 inches
Peso del generador	77kg / 169,4lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

