

NG3000A Pureza ultra alta Generador de nitrógeno

Aplicación: Soluciones especializadas


Número de pieza: 10-1210 (230v) - UK | 10-1220 (230v) - EU | 10-1230 (110v) - US



Description

El NG3000A entrega de nitrógeno de ultra alta pureza en los límites de detección de trazas a través de la tecnología de adsorción por cambio de presión en los flujos de 3 l / min para aplicaciones GC .

- **Tipos de gas** Nitrógeno
- **Compresor integrado:** Si

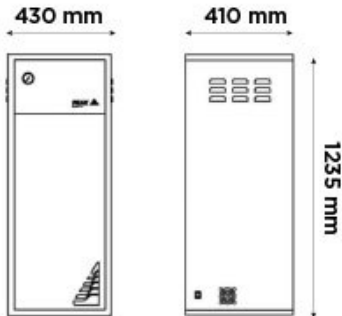

Key Features

- Adecuado para aplicaciones como GC, DSC, DMA
- Adecuado para aplicaciones que requieren una pureza ultra alta
- Disponible con o sin compresor externo
- Contiene una columna de tamiz molecular de carbón autorregenerante
- El gas se suministra a pedido para que el generador funcione según su horario
- Garantía integral in situ de 12 meses

Technical Specifications

NG3000A Pureza ultra alta Generador de nitrógeno	
Tipo de gas	Nitrógeno
Flujo de gas máximo:	3
Presión máxima de salida:	80psi/5,5bar
Presión máxima de salida:	80psi/5,5bar
Pureza máxima:	99.9995%
Voltaje:	115 / 230 ± 10% VAC
Frecuencia:	50 / 60 Hz
Corriente:	3.6 Amps / 8 Amps
Acreditaciones:	CE
Tamaño (HxAxD) en mm	1235 x 449 x 406,5 mm
Tamaño (HxAxD) en pulgadas	48,6 x 17,7 x 16 inches
Peso del generador	88kg / 193,6lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

